

Combination of interdental brush with toothpick!

This is to take out the food particles lodged between your teeth or to remove plaque in the space between the teeth that cannot be removed by brushing.

It removes plaque and bacterial membranes left between your teeth after brushing.

This is recommended for people with wide gaps in the root area of the teeth or for people wearing braces.



Specifications


Product weight : 22g Container weight : 7g
 Lid weight : 1g One piece weight : 0.2g
 Product length : 65mm
 Product thickness : 0.6t (toothpick part) 0.45t (interdental brush)
 Material of the container / Cover : Polypropylene (P.P.)
 Material of the product : Polyoxymethylene (P.O.M)

Other details


1. Food particles that have not been removed can be easily removed after brushing.
2. Put the interdental brush part between the teeth and move it back and forth to remove any food particles.
3. Remove any residues with a toothpick on the other side.
4. You can put a interdental brush between your teeth and move it back and forth to remove a toothbrush that cannot be removed.
5. If the gums bleed while using it, stop using it immediately.
6. If the bleeding continues, consult dentist.
7. Due to the plastic properties, do not heat the product directly or store it in a firearm or hot place.

Combination of interdental brush with toothpick!

Usage

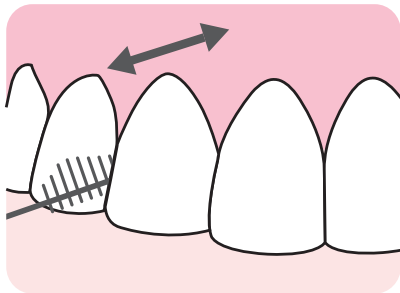

Place the interdental brushin between the gums and move it back and forth to clean the food particles.

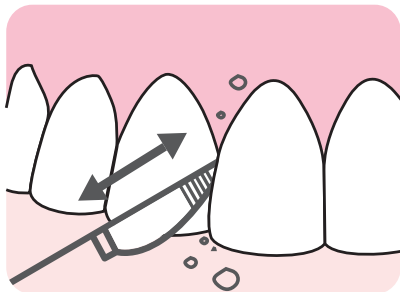

Remove food particles and tartar with saw blade-like parts.

Notice

If a person with narrow teeth inserts this too much between them, your teeth or gums can get hurt. If gum disease or bleeding continues during use, consult dentist.

